### Vetus 8

Veterinary Diagnostic Ultrasonic System

# Excellent performance for your trust










P/N:ENG-Vetus 8-210285X8P-20200121 ©2020Shenzhen Mindray Bio-Medical Electronics Co.,Ltd. All rights reserved.





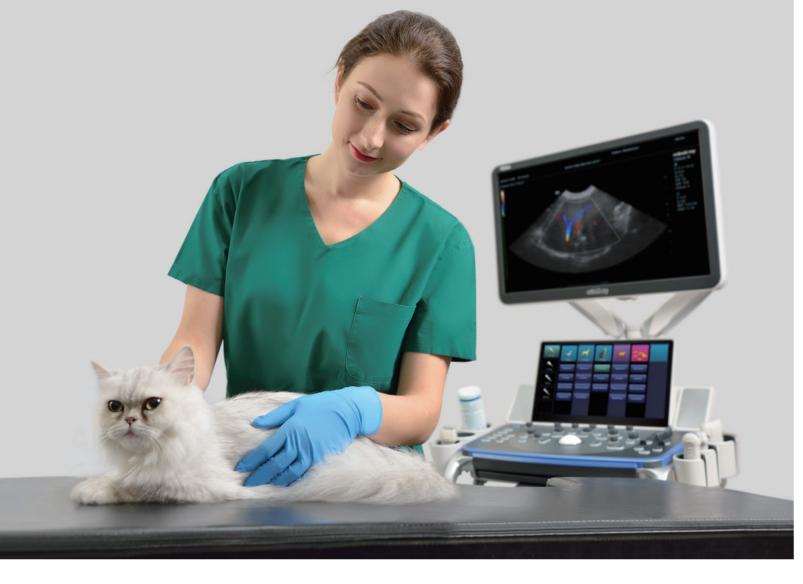

# Vetus 8

Mindray Vetus Series offer professional veterinary solutions catering to a wide variety of species ranging from small pets to large farm animals.

Vetus 8, the premium veterinary ultrasound imaging system, adopts the most cutting-edge ultrasound technology to provide clear and precise image according to different animal species.

Moreover, packaged with the simple workflow, dedicated design and complete veterinary solution make it the ideal ultrasound system for your everyday practice.

### **Clear & Precise Images**

- Professional veterinary image presets according to different animals' body characteristics
- Animal species are subdivided to different categories by weight and body size for precise diagnosis
- 23.8 inch High Definition LED monitor provides clear tissue details

### Simple & Smart Workflow

- Transducers and the according exam modes are displayed in the same interface, select or switch within 1 second
- Recent 4 modes are presented on screen, one-key to switch the frequently-used exam mode
- 13.3 inch touch screen supports powerful multi-touch gesture operation




### **Powerful Intelligent Tools**

- Auto Volume Control System helps to prevent frightening the pet animals
- Auto measure tools: Auto PW, Auto EF maximize the productivity for easy-accurate-scanning
- Automated workflow guide iWorks standardize and simplify the Abdomen/Cardiology workflow, let you be more focus on the animal diagnosis

### iTouch

iTouch is a one-button auto image optimization solution, which is able to optimize the imaging quality automatically, including in B/PW mode.

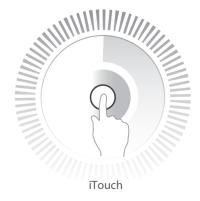





Liver Feline




Long Axis View Cardiac Canine

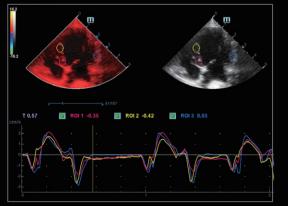


Cardiac TDI QA Canine

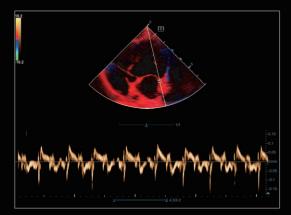


Cardiac TVD Canine

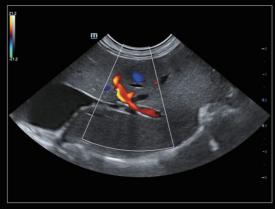


Hepatic Flow Feline

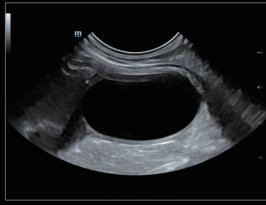

Bladder Feline

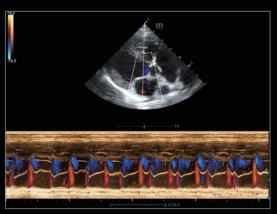


Cardiac Color M mode Canine

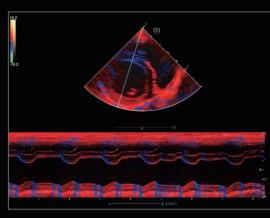

Cardiac TVM Canine

## Ergonomics dedicated for veterinary practice

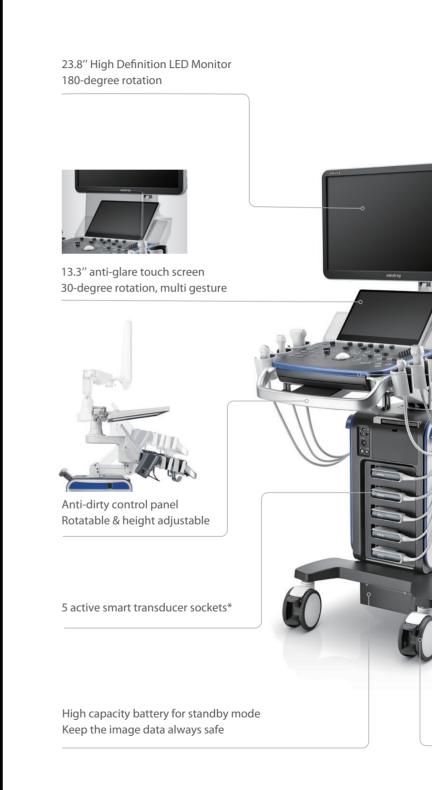

\*4 sockets is standard configuration, 5 sockets is purchasable optional.



Dual-Wing Floating Arm Move anywhere you want



Sliding keyboard with designed silicon cover

Designed for cable management Help keep cables off the ground

Stable wheels with small footprint Anti hair-wrapped, 360° rotation

### **Small Animal Solution**

Cat, small dog and exotic small animals

#### **Abdomen & Reproduction**



#### C11-3s

Classical Micro convex array transducer for abdomen

- Dedicated design: streamlined and smart size for small animal abdomen
- One probe solution: Support both abdominal and basic heart ultrasound
- Excellent performance: High penetration with wide field of view (Extended FOV: 140°)

#### Superficial & Musculoskeletal



#### L12-4s

High frequency linear array transducer for small parts

- iNeedle: Biopsy needle enhancement and guidance
- iScape View: Real-time Panoramic Imaging to extend the field of view
- HR Flow: High resolution flow to provide more details of animal tissue
- Elastography imaging: displaying the tissue stiffness for precise evaluation

#### Cardiology



P8-2s Mid/high-frequency phased array for mid and small-sized species

P10-4s high-frequency phased array for small-sized species

- Free Xros M-Anatomical M module: up to 3 sample lines
- TDI-Tissue Doppler Imaging with complete 4 modes (TVI/TVD/TVM/TEI): wall motion analysis for myocardial function
- TDI QA-Tissue Doppler Imaging with Quantitative Analysis: analysis of myocardium motion with strain/strain rate
- Auto EF: Automatic animal heart Ejection Fraction computation
- Stress echo: Evaluate segmental myocardium movement with increased workload
- TTQA-Tissue Tracking with Quantitative Analysis: track myocardial contraction by detection of the 2D speckle patterns, provide LV regional movement abnormities diagnostic

### Large Animal Solution

big size dog, horse, cattle, sheep and exotic big animals

#### Reproduction

#### 6LE5Vs

Intrarectal linear array transducer for equine and farm animal reproduction exam

- 3 m cable tailed for large animal scenario
- Anti-liquid, anti-dirty and easy-clean design
- High resolution with wide field of view designed for IVF

#### Abdomen



#### C6-2s

Convex array transducer for abdominal and reproduction scanning in big animals

- iNeedle: Needle enhancement and guidance for biopsy
- iScape View: Real-time Panoramic Imaging to extends the field of view

#### Cardiology

#### P4-2s

low-frequency phased array transducer for medium-sized & large species

- Free Xros M-Anatomical M module: up to 3 sample lines
- on the image
- analysis for myocardial function

\*More transducers can be supported, please contact your local sales representative for further details.

• Free Xros CM-Curved Anatomical M module: the sample line manually depicted anywhere

• TDI-Tissue Doppler Imaging with complete 4 modes (TVI/TVD/TVM/TEI): wall motion